

**Alc
Smart**

Breath Alcohol Tester
AT6000 **Operation Manual**



Hanwei Electronics Corp.,Ltd.

169 Xuesong Road, National Hi&New Technology Zone, Zhengzhou China 450001

Tel: 86-371-67169070 67169080

Fax: 86-371 67169090

Email: hwsensor@163.com

<http://www.hwsensor.com>

HW14112504WS

Preface

Reminding Before your testing.

Alco Smart Breath tester should be used only to give an indication of the possible presence of alcohol in the breath/blood. You should not rely upon it as the sole basis to determine intoxication or whether it is safe to drive a vehicle, operate equipment, or engage in dangerous activities.

Everyone has different body responds to alcohol consumption and his testing result only as reference, not a subject standard for consequent decision.

The manufacture, importer or distributor takes no responsibility whatsoever for the use of this product for any reason. This product must not be used as a tool for determining whether a person is able to operate a motor vehicle or device legally or safely. The intake of any alcohol will impair reflexes and judgment to operate motor vehicle.

Before your testing, please view the manual carefully and follow the instruction strictly.

General Introduction

AT6000 breath alcohol tester is kind of tester designed to measure concentration result of the breathed alcohol in the human body. This device adopts advanced flat surfaced alcohol sensor which has excellent sensitivity and reproducibility, fashionable and portable design make it more convenient for personal use. When the alcohol content exceeds the preset limited level, this device would send an audio and video warning to remind your safety.

Main Feature:

Advanced flat surfaced alcohol sensor

Quick response

SMD assembling, stable performance

Smart MCU control

Direct testing process LCD indication

Digital LCD display with light blue backup

Portable and fashion design

Audio warning beyond pre-set limit

Sensor fault self checking

Battery saved design, low voltage indication

Technical Data

Sensor type: Flat surfaced alcohol sensor

Detection Range: 0.00 – 0.20%BAC; 0.00~2.00g/L;

0.00 – 2.00‰BAC; 0.00 -1.00 mg/L; 0.00 – 200mg/100ml

Alarming Level: 0.05%BAC; 0.50g/L; 0.50 ‰BAC; 0.25mg/L; 50mg/100mL

Accuracy: $\pm 0.02\%$ BAC at 0.05%BAC

Working Voltage: DC4.5V (3xAAA) Batteries; Working Current: $\leq 120\text{mA}$

Working Environment: Temperature $0^{\circ}\text{C}\sim 50^{\circ}\text{C}$

Relative Humidity $\leq 95\%$ No Dews

Display: 3 digits LCD display with light blue backup
 Dimension: 103×65×27mm (L×W×H), ≤64 g Battery excluded
 Battery expectancy time: ≥200Times


1. Structure and Function Guide

1.1 Structure Guide



1	Mouthpiece
2	Buzzer hole
3	Switch Button
4	Input for external power
5	LCD Screen
6	Vent hole

1.2 Function guide

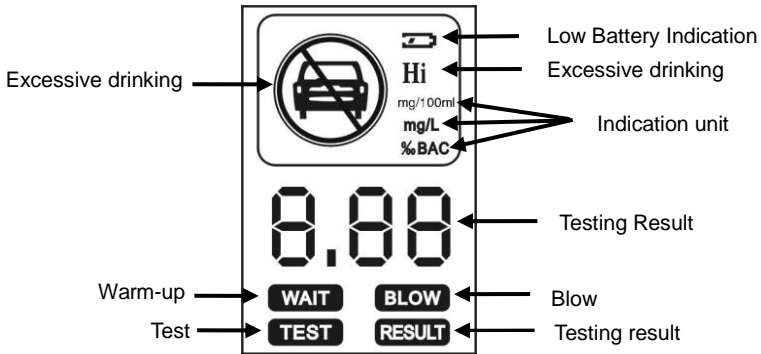

Picture 2

8.8 8	Alcohol content
"C" shows circularly on the LCD	ready for blow
"—" shows circularly on the LCD	Analyze the concentration, please wait
FFF	sensor fault

Digital Indication:

2. Operation Instruction

- 2.1. Slide down the battery compartment on the back side, insert 3pcs AAA 1.5V batteries inside according to battery polarity indication and close it.
- 2.2. Press the switch button and last for 1s, the tester would be power on with a buzzer brief

ringing, 'Wait' shows on down side of screen, meantime it shows 100 to 000 countdown, please wait now, otherwise its accuracy will be not good.

Until "Wait" disappears, tester will send out a ring, meantime shows "BLOW" and "C" gliding, now the test could begin.

Remarks: Its countdown speed will depend on concentration of last test, if concentration is high, it will be slow. If concentration is low, the countdown will be very quick.

2.3. Have a deep breath before testing, and blow directly to gas entry when "C" circularly shows, stop blow when heard its buzzer ring.

2.4. After blowing finished, the LCD screen will show "●" circularly, meantime shows "TEST", this means the tester is analyzing the result, do not blow at this time. After analyze the screen will show the testing results.

It will turn off automatically after keep results for several seconds.

2.5. If a new testing needed, press the switch button again and repeat the above step 4.2~4.4.

3. Notification

3.1. Avoid any fall or strong shock.

3.2. If noise gas with high concentration existed, may the tester won't work normally.

3.3. If testing is done under low voltage, certain error will be existed between the real value.

3.4. To ensure the testing result, please wait 15 mins to take the testing after your drinking.

3.5. The result may be unreliable for your first test if long time storage, please try second time.

3.6. Do not keep the tester in the environment of corrosive gas (Chlorine etc) for use or deposit, also in other bad surroundings.

3.7. Do not blow smoke to tester directly, especially the smokers.

3.8. When in consecutive testing, if blow immediately after drink, this concentration is very high, there will be some rudimental in mouthpiece. Please wait for at least 1 min to take second test, or, blow clean air for 3 seconds to clean out its mouthpiece.

3.9. After Long time use, there may be dirt on the tester, please use clean cloth to wipe off. Do not use any hard object or the solvent with any Corrosive ingredient.

Remarks: Common Fault and Solving solutions

Fault	Possible reason	Solution
NO display on LCD	Incorrect battery installation	Insert battery correctly according to polarity
	Low battery seriously	New battery replacement
	Circuit fault	Contact distributor
No response to detection gas	not complete Warm-up	Waiting complete Warm-up
	Circuit fault	Contact distributor
FFF display then power off	Sensor fault	Contact distributor
Low voltage display then power off	Low battery seriously	New battery replacement

To keep the continuous improvement of the products, we reserve the right to improve design without prior notice.